



STARJET *

■ stop limitations

A typical scenario - you complete the supragingival air-polishing with your turbine-adapted handheld unit and your next step is subgingival cleaning but your unit is not intended for this use.

► start flexibility

With the STARJET, these limitations no longer exist. With the flexibility of choosing the appropriate powder, you can easily switch from supra- to subgingival by simply turning the selection ring between prophy and perio. Not flexible enough? You can also select your preferred spray head, with three options.



120°



PERIO



90° (optional)




► 3 different spray heads

Always best access - 90° and 120° angled spray heads for supra- and subgingival use in periodontal pockets up to 5 mm depth.

The spray head PERIO equipped with the disposable sterile Subgingival Perio Tip provides optimal access in pockets deeper than 5 mm.

Just one click connects the spray head safely to the STARJET, thanks to a dedicated security connection system.

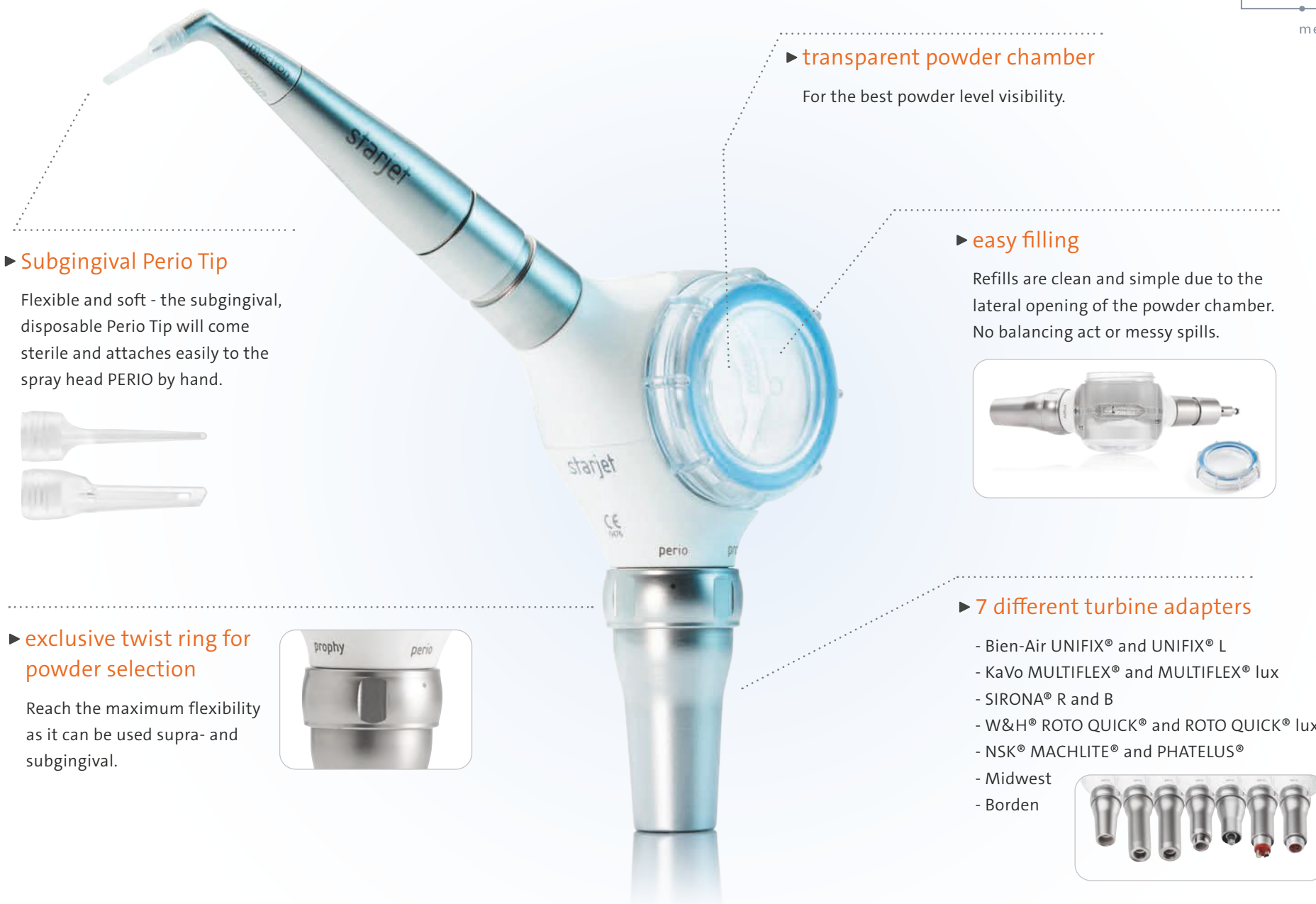

► Subgingival Perio Tip

Flexible and soft - the subgingival, disposable Perio Tip will come sterile and attaches easily to the spray head PERIO by hand.



► exclusive twist ring for powder selection

Reach the maximum flexibility as it can be used supra- and subgingival.



► transparent powder chamber

For the best powder level visibility.

► easy filling

Refills are clean and simple due to the lateral opening of the powder chamber. No balancing act or messy spills.



► 7 different turbine adapters

- Bien-Air UNIFIX® and UNIFIX® L
- KaVo MULTIFLEX® and MULTIFLEX® lux
- SIRONA® R and B
- W&H® ROTO QUICK® and ROTO QUICK® lux
- NSK® MACHLITE® and PHATELUS®
- Midwest
- Borden





AIR-POLISHING POWDERS

■ stop too much choice

The number of powders for dental use is countless - how on earth could you possibly find exactly the one you need to achieve the best results?

► start always the right choice

Here's a tip: just go for the well-known Mectron quality. Our specially developed air-polishing powders not only come with a hydrophobic effect, which prevents clogging perfectly – but, more importantly, they also guarantee perfect air-polishing results.⁹⁻¹⁶



► PROPHYLAXIS POWDER SOFT⁺

gently and effectively removes supra-gingival biofilm and stains on the enamel surface (sodium bicarbonate-based - average particle size 40 µm)

Indications:

- Removal of stains and extrinsic discolorations due to smoke, wine, tea and coloring agents
- Cleaning of orthodontic brackets















► **PROPHYLAXIS POWDER
SENSITIVE⁺ – GLYCINE POWDER**

gently removes supra and subgingival
biofilm (glycine-based - average
particle size 25 µm)

Indications:

- Mild pigmentation and stain removal
- Periodontal pocket detoxification
- Periodontal maintenance
- Implant maintenance and treatment
of peri-implantitis
- Orthodontic patient maintenance
- Gentle cleaning of restoration materials
and comfortable treatment on exposed
root surfaces
- Reduction of bacteria colony-forming
units (CFU)
- Can be successfully used for the
treatment of each type of patient

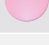INDICATIONS	SOFT ⁺	SENSITIVE ⁺
SUPRAGINGIVAL BIOFILM REMOVAL		
REMOVAL OF STAINS AND PIGMENTATIONS		
SUBGINGIVAL BIOFILM REMOVAL		
PERIODONTAL POCKET DETOXIFICATION		
IMPLANT MAINTENANCE		
PERI-IMPLANTITIS TREATMENT		
FIXED ORTHODONTIC APPLIANCES		
CLEANING OF PROSTHESIS AND RESTORATION MATERIALS		
RESTORATION MAINTENANCE	