



27 CLINICAL EVALUATORS

468 TOTAL USES

86% CLINICAL RATING

Key features: Self-curing composite • Gap free • Unlimited depth of cure
• Two steps • Available in syringe or capsule delivery

Description

Stela properties:

- Self-curing composite
- Unlimited depth of cure
- Contains fluoride, calcium and strontium

Photo by Dr. Matthew Miller



Core build-up using Stela composite

Indications

- Class I, II, III, & V
- Core Build Ups
- Base/Liner
- Sealer for endodontic access openings

What procedures did you use this product for?

Answer Choices	Responses
Class I restoration	52%
Class II restoration	52%
Class III restoration	7%
Class V restoration	30%
Core build up	82%
Base/liner	33%
Sealing endodontic access openings	41%



Clinical Tips

- Use proper isolation as this is a self-cured material.
- For Class II restorations, use a metal matrix so you can sculpt the material.

Unique Attributes

- Self-cured primer and composite
- Only two steps: prime then restore
- Gap free: primer initiates curing process at restoration interface
- Unlimited depth of cure
- 15 second preparation, sets in 4 minutes
- High flexural and compressive strength
- Can be used on all tooth shades in non-esthetic zones
- Anterior for A2/A3 Patients
- BPA and HEMA Free

"FAST SETTING
TIME, GOOD
ESTHETICS."

Evaluators' Comments

"The self-curing capability is great, especially in deep access cavitations where lack of confidence is found in light curing."

"It did not take more time than incremental additions of light-cured composite."

"Large depth of cure, plenty of time to shape and sculpt, minimal steps."

"Very easy to polish. Material wasn't chalky like other materials."

"Set time was inconsistent and took too long for my taste."

Consultants who would:

63% Recommend to a colleague


Consultants who would want to stock in office:

11% Yes, instead of current product

41% Yes, in addition to current product

18% I might want to order this product for certain cases

Evaluation Summary:

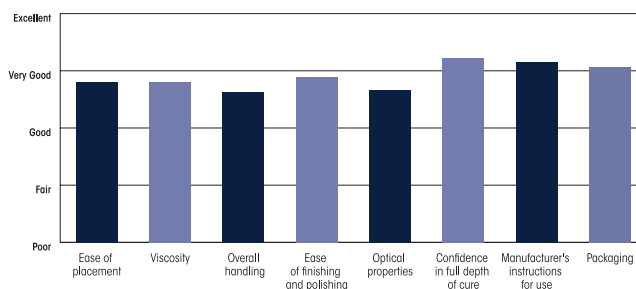

Compared to Competitive Products:

